DISPO PREP



DISINFECTANT CAP

Smart Protection. Continuous Defense.

Your First Line of Defense Against Catheter-Related Infections

SURGEINE HEALTHCARE


DISINFECTANT CAP

Disinfectant Cap - Passive Protection for Active Safety

Surgeine's disinfectant caps are single-use, alcohol-impregnated devices designed to provide continuous passive disinfection of needleless IV access ports. Each cap contains 70% isopropyl alcohol (IPA), ensuring effective microbial reduction between line accesses. Brightly colored for visual confirmation, these caps support standardised catheter maintenance protocols and play a critical role in reducing central line-associated bloodstream infections (CLABSIs). Easy to apply and universally compatible, they are a key component in any infection prevention strategy.

Product Features:

- **Continuous Disinfection:** Maintains antiseptic barrier for up to 7 days if left in place
- **70% IPA Foam Core:** Effectively saturates the septum and connector threads
- **Color-Coded Cap:** Enhances visibility and protocol compliance
- **Universal Fit:** Compatible with all ISO 80369-7 standard luer-access needle-free connectors
- **Sterile and Latex-Free:** Safe for all patient populations

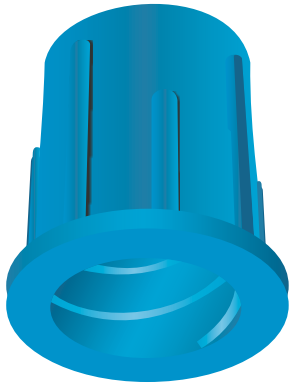

Key Benefits

- Significantly reduces CLABSI and CRBSI risk
- Provides passive disinfection for up to 7 days
- Safe for use across central lines, PICCs, and peripheral Ivs
- Offers a visual indicator for cap presence and protocol adherence
- Ensures compliance with standardized infection prevention practices

Why Choose Us

- *First in India to manufacture in-house*
- *Supplied to over 12 countries across Europe, the Gulf, Africa, and ASEAN*
- *Fully sterile production validated with ISO 13485 & MDSAP standards*
- *Optional experience center for evaluation and live product demonstration*

SURGEINE HEALTHCARE

Clinical Endorsements

The Infusion Therapy Standards of Practice (INS) recommends the use of passive disinfection caps with 70% IPA due to their association with significantly reduced rates of catheter-related bloodstream infections..

The Society for Healthcare Epidemiology of America (SHEA) supports disinfectant caps as a proven tool to enhance infection control protocols across care settings.

Regulatory Compliance & Standards

- *Compliant with ISO 80369-7 standards for connector compatibility*
- *Aligned with INS and SHEA guidelines for infection prevention*
- *Clinically validated for maintaining disinfection barrier up to 7 days*

Image-1

PEEL

Remove protective cover

Image-2

PUSH

Twist: Attach cap clockwise

Image-3

PROTECT

Leave in place until next access

Image-4

PURGE

Twist counterclockwise and discard

Packaging Formats

Code	Form	Qty/Box
3051	Single Pieces	250 Caps
3052	Strips of 5	500 Caps
3053	Strips of 10	250 Caps



Get in touch with us:

SURGEINE HEALTHCARE (INDIA) PRIVATE LIMITED

A-310, Urban Square, Sector-62, Gurugram, Haryana - 122 002, INDIA

+91 91118 89 777 **sales@surgeine.com** **www.surgeine.com**

DISCLAIMER: The information supplied in this document is for guidance only and should not be construed as a warranty. All implied warranties are expressly disclaimed, including without limitation any warranty of merchantability or fitness for use. All users of the material are responsible for ensuring that is suitable for their needs, environment and end use. All data is subject to change as **SURGEINE HEALTHCARE (INDIA) PRIVATE LIMITED** deems appropriate.

Scan QR Code on Your Smart Phone or Ipad to Save A Digital Version